

ENTityXL LED NasoView Nasopharyngoscope

Material compatible with STERIS Revital-Ox™ Resert®, Rapicide OPA/28, Metricide OPA Plus



The ENTityXL Nasopharyngoscope is now material compatible with more High Level Disinfectants

The ENTityXL LED NasoView Nasopharyngoscope frees you to perform nasal endoscopy and laryngoscopy virtually anywhere. Our patented, integrated LED light source eliminates cumbersome external light sources while providing exceptional visibility.

The ENTityXL has been tested in Optim's testing lab for material compatibility with the STERIS Revital-Ox™ Resert® HLD Disinfectant system, Rapicide OPA/28 and Metricide OPA PLUS for 1,000 cycles. This is in addition to passing the STERIS Device Testing Program's rigorous material compatibility evaluation for the STERIS Resert Disinfectant system.

Specifications

Flexible Length	30cm (11.8")
Overall Length	53cm (20.9")
Shaft Diameter	3.6mm
Weight	11.6 oz
Articulation Range	135 degrees up / down
Depth of Field	5mm to 50mm
Field of View	70 degrees
Battery	Lithium Ion

Compact, portable, all-in-one design


- No external light source or lamp replacement
- Hybrid aluminum handle
- Fully sealed and completely immersible
- Bring the scope to the patient
- Use in the office / bedside / hospital

Better economics

- See more patients per day, increase revenue
- Reduced capital expenses and ongoing maintenance
- No lamp or light guide replacements

Superior lighting

- Patented LED light coupling
- Enables clear visibility of the true vocal chords
- Variable light intensity

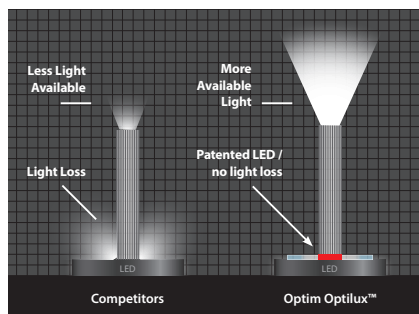

Optilux™

Optim's unique LED light source to fiber light guide coupling technology, Optilux™, in the EntityXL NasoView brings unmatched illumination to the end of the scope, where it counts.



64 Technology Park Road, Sturbridge, MA 01566
Toll-free: 800.225.7486 Worldwide: 1.508.347.5100
Fax: 508.347.2380
Email: sales@optim-llc.com
Website: www.optim-llc.com



Optim is ISO 9001:2008 / ISO13485:2003 Certified